
# Your PDF Guides

You can read the recommendations in the user guide, the technical guide or the installation guide for BRAUN IPL VENUS NAKED SKIN BY BRAUN. You'll find the answers to all your questions on the BRAUN IPL VENUS NAKED SKIN BY BRAUN in the user manual (information, specifications, safety advice, size, accessories, etc.). Detailed instructions for use are in the User's Guide.

**User manual BRAUN IPL VENUS NAKED SKIN BY BRAUN**  
**User guide BRAUN IPL VENUS NAKED SKIN BY BRAUN**  
**Operating instructions BRAUN IPL VENUS NAKED SKIN BY BRAUN**  
**Instructions for use BRAUN IPL VENUS NAKED SKIN BY BRAUN**  
**Instruction manual BRAUN IPL VENUS NAKED SKIN BY BRAUN**

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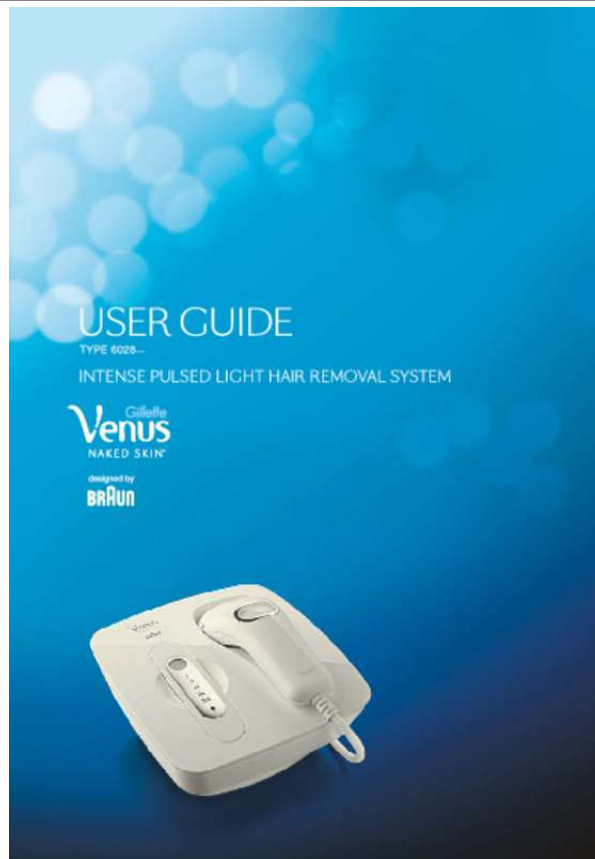If you are unsure if the device is appropriate for you to use, please call the toll-free number located on page 41 of this user guide or consult your physician.

WARNING – DO NOT USE: If the skin to be treated is darker than skin tone (III) shown on the Skin Tone Chart.

[You're reading an excerpt. Click here to read official BRAUN IPL VENUS NAKED SKIN BY BRAUN user guide](http://yourpdfguides.com/dref/5484263)  
<http://yourpdfguides.com/dref/5484263>

Darker skin tones absorb too much light energy during treatment, which may injure (e.g. burn, discolor or scar) your skin. As the face can be more sensitive, use the device on other body areas before treating upper lip area, chin or jaw. For specific facial usage instructions, see chapter "Facial Treatment".  
Welcome to Gillette Venus Naked Skin® Welcome to a premium beauty experience with your new Venus Naked Skin® Intense Pulsed Light Hair Removal System. With Venus Naked Skin®, you can achieve smooth skin every day, so that every season can be bikini season. So get started... you'll love the look and feel of your continuously smooth skin.

And that's a beautiful thing! III III DO NOT USE 22 For more information and to download the free Naked Skin® mobile application, please visit [www.braun.com](http://www.braun.com) 23 Contents Who Should Not Use Gillette Venus Naked Skin® (Contraindications).....

.....

.....

.. 24 Before You Start (Warnings and Precautions).....

.....

.....

.....

.....

.....

.....

.....

. 24 Warnings.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.. 24 Precautions...

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

..... 25 Electrical and fire safety.

.....

.....

.....

.....

.....

.....

.....

.....  
.....  
.....  
.....

.....  
*.. 26 Possible Side Effects...*

.....  
.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

*27 Possible side effects include:.....*

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....

*27 What to do if you have a side effect....*

.....  
.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

*27 Introduction to Gillette Venus Naked Skin® .....*

.....  
.....  
.....

.....  
.....  
.....  
.....

*.. 28 How It Works.....*

.....  
.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

..... 29 Check Your Skin Tone.

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

..... 29 Treatment.....

.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

31 Facial Treatment (women only).....

.....  
.....  
.....

.....  
.....  
.....  
.....

.....

.....  
.....  
.....  
... 33 Patch Test.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

. 35 Treatment Schedule.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.. 36 Start-up Phase.....

.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....

.... 36 Maintenance Phase.

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

.....  
.....  
.....  
.....  
.....  
..... 36 *Troubleshooting.*

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

.. 37 *Technical Specifications..*

.....  
.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....

. 38 *Customer Support.....*

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....

..... 39 *Warranty.....*

.....  
.....  
.....

.....  
.....  
.....  
.....

.....  
.....  
.....  
.....  
.....  
.....  
.....  
.....

..... 39 25 Who Should Not Use Gillette Venus Naked Skin® (Contraindications) • **DO NOT USE** the device on scalp, eyebrows or anywhere around the eye or eye D socket. Use in these areas may permanently damage eyes. • **DO NOT USE** if you have dark skin tones (see skin tone chart on page 22). Dark skin D may absorb too much light energy and injure (e.g. burn, discolor or scar) your skin. • **DO NOT USE** if you have a history of skin cancer or precancerous lesions (e.g., D nevi or a large number of moles). Before you start, consult your doctor if you take medication for any chronic medical condition.

• **DO NOT USE** if you are pregnant, lactating, or under the age of 18 as this device has D not been tested with these individuals. – **DO NOT** a tattoo or permanent makeup On on the area you wish to treat. – **DO NOT** you have a known sensitivity to If sunlight (photosensitivity) or are taking medication that makes the skin more sensitive, e.g., Retin A, Accutane and/or other topical retinoids.

– **DO NOT** if you have had a skin peel I treatment on the area you wish to treat. • **DO NOT** if you are taking any kind of medication I regularly and over a long period of time, please consult your physician for any potential impact on skin sensitivity to light. • **DO NOT** his device is not intended for use T by persons under 18 years of age or persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. • **DO NOT** children should be supervised to C ensure that they do not play with the device. When not in use, always unplug the device. The device should be stored securely out of the reach of children. • **DO NOT** replace the gel with other gels, D creams or lotions. Use only the gel specifically designed for use with the Venus Naked Skin® device. • **DO NOT** void contact of gel with the eyes. In A case of accidental contact between the gel and eyes, rinse with water.

If skin irritation or rash develops, discontinue use of the gel and the device. If symptoms persist, consult your physician. • **DO NOT** regularly check the device and cords R for visible signs of damage. In case of damage or cracks, stop using the device and have it checked and repaired by an authorized service repairman. Call the toll-free number on page 41 for assistance. Unqualified repair work can lead to extreme hazards for the user. • **DO NOT** void any activity or skin products that A may irritate your skin for 24 hours after treatment. This includes: – Using hot tubs or saunas **DO NOT** – Wearing tight-fitting clothing **DO NOT** – Using bleaching creams, perfumed **DO NOT** products or peeling products 27 Before You Start (Warnings and Precautions) WARNINGS • **DO NOT USE** Venus Naked Skin® if **DO NOT** your skin is darker than the dark skin tone (IIII) shown on the skin tone chart or if the setting (IIII) red light on the Skin Tone Sensor is lit. Because dark skin absorbs more light energy, treating dark skin with Venus Naked Skin® may cause discomfort/pain or may injure (e.g. burn, discolor or scar) your skin. • **DO NOT USE** Venus Naked Skin® **DO NOT** if your skin has been artificially or naturally tanned recently. Because tanned skin is darker, it is at a higher risk of damage when using Venus Naked Skin®. Treating tanned skin with Venus Naked Skin® may cause discomfort/pain or may injure (e.g. burn, discolor or scar) your skin. Avoid unprotected exposure to direct sunlight after your Venus Naked Skin® treatment. Treatment with Venus Naked Skin® may make your skin more sensitive to the sun, resulting in side effects such as burning and blistering. Use sunscreen (SPF 15 or greater) or cover the skin. • **DO NOT USE** Venus Naked Skin® on: – Red hair – Grey or white hair – Very blonde hair Venus Naked Skin® is not effective on these hair colors. You will be exposed to the risks of this product but not the benefits. • **DO NOT USE** Venus Naked Skin® D directly on nipples, genitals or around the anus. These areas may have a darker skin color and/or greater hair density and using the device in these areas may cause discomfort/pain or injure (burn, discolor or scar) your skin. • **DO NOT USE** Venus Naked Skin®: – **DO NOT** n dark brown or black spots such O as freckles, birthmarks, moles or warts in the area you wish to treat. – **DO NOT** if you have chronic skin disease I (e.g., psoriasis or vitiligo) or damaged skin (e.g., sunburn, cuts, open wounds or active infections) in the areas you wish to treat. – **DO NOT** if you have varicose veins in the I area you wish to treat.

If you have any of the above, Venus Naked Skin® may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin color changes (hyper- or hypo-pigmentation), or scarring. • **DO NOT** look directly into the glass D filter on the handset or attempt to activate the handset toward the eyes. To protect the safety of your eyes, and those of individuals nearby, the Venus Naked Skin® handset will only activate when both contact bars, on the bottom of the handset, are placed against the skin. If at any time the handset activates when both bars are NOT in contact with the skin, STOP use of the device immediately and contact your retailer. • **DO NOT USE** Venus Naked Skin® more D than once a week on the same area. It will not speed up the results and may increase the chances of injury to your skin.

[You're reading an excerpt. Click here to read official BRAUN IPL VENUS NAKED SKIN BY BRAUN user guide](http://yourpdfguides.com/dref/5484263)  
<http://yourpdfguides.com/dref/5484263>

**PRECAUTIONS** • Do not use tanning including outdoor sun A bathing, tanning treatments in a salon before and during treatment, or use of a tanning lotion or cream before and during treatment as the device may not work for you. Further, darker skin may absorb too much light, which may cause discomfort/pain or may injure (e.g.

burn, discolor or scar) your skin. 26 The above may cause irritation, such as itching and redness, in the treated areas. **ELECTRICAL AND FIRE SAFETY** As with any electrical device, certain precautions must be taken in order to ensure your safety. **DANGER:** Never open the Venus Naked Skin® device, as this may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury. • Do not use if device is damaged; for example, cracked base, handset cable damage (internal wires visible), cracked or broken glass on the handset, etc.

Internal parts contain electric charges that are hazardous. Electric shock may result from using a damaged device. • Do not use Venus Naked Skin® near D bathtubs, showers, basins or other vessels containing water. If Venus Naked Skin® becomes wet, it should not be used. Electric shock may result from using a wet device.

• Do not use if either the base unit or D handset becomes uncomfortably hot to touch. The device may overheat and pose an electric shock risk or fire hazard. • Do not use a power supply other than provided Cyden Type 6027... supply for use with Venus Naked Skin® Type 6028... This may damage the unit and pose an electric shock risk or fire hazard • Do not use petroleum-based or D flammable cleaning agents as these pose a fire risk. • Do not use abrasive cleaning D agents. These may damage the device or scratch the glass window on the handset. Damaged devices must not be used. Electric shock may result from using a damaged device. • Do not block the cooling vents on D the handset. The device may overheat and pose an electric shock risk or fire hazard.

**NOTICE:** • Venus Naked Skin® is designed for use V where the main supply voltage is 100 V to 240 V AC, 50/60 Hz. Check that the actual supply voltage is within this range. • Always unplug Venus Naked Skin® A when not in use. **POSSIBLE SIDE EFFECTS** Venus Naked Skin® is a powerful device. Some people will experience slight discomfort (e.g., heat or redness of the skin) when using Venus Naked Skin® – this is normal. You may feel heat and tingling when activating the light pulse, followed by redness in the skin. This type of reaction usually disappears within 24 hours. Adverse reactions are usually immediate or appear within 24 hours.

In rare cases, they may take up to 72 hours to appear. **POSSIBLE SIDE EFFECTS INCLUDE:** • Discomfort/pain in the area being treated D • Itching of the treated area I • Heat and redness of the skin that H continues after the first 24 hours In rare cases: • Itching or burning of the skin B • Swelling and / or bruising S • Changes to the pigmentation (color) of C the skin in the area treated **WHAT TO DO IF YOU HAVE A SIDE EFFECT** Side Effect Warm feeling or tingling sensations during treatment which typically disappear after a few seconds to a minute and decrease with continued use. Discomfort or pain is intense during treatment or persists after a treatment. Skin redness during or after treatment which disappears within minutes to several hours. Skin redness which does not go away within 24–48 hours after treatment.

**How to Assess and React** This is expected and is normal for all IPL (Intense Pulsed Light) hair removal device treatments. You can keep on using the device as instructed. Stop using the device and consult your physician before using it again. This is expected and is normal for all IPL (Intense Pulsed Light) hair removal device treatments. You can keep on using the device as instructed once skin redness has disappeared.

Stop using the device and consult your physician before using it again. In Very Rare Cases: Swelling and redness around the hair follicle that disappears within two or three days. Swelling and redness around the hair follicle that does not disappear within two to three days. Temporary changes in skin color (lightening or darkening). Increase in hair growth in areas of treatment (more likely to be seen with subjects of Mediterranean descent). You can keep on using the device as instructed once swelling and redness around the hair follicle has disappeared. Stop using the device and consult your physician before using it again. If your skin color changes, stop using the device and consult your physician. If this increase persists over several treatments, stop using the device and consult your physician. 28 29 Introduction to Gillette Venus Naked Skin® The Venus Naked Skin® device is powered by Intense Pulsed Light (IPL), derived from professional light technology that is used by dermatologists/professionals worldwide for the treatment of unwanted hair.

The Venus Naked Skin® System includes: 1. Base unit 2. Handset (connected to base unit) 3. Detachable Skin Tone Sensor (see page 31 for details on the Skin Tone Sensor) 4. Activation button 5. Power supply 6. Treatment window in the front of handset 7. <sup>TM</sup> What to Expect You could experience noticeable hair reduction after 6 treatments and long lasting hair removal after 6 – 12 treatments. Continue to treat once every two months to maintain results and to enjoy continuously smooth, beautiful skin. 7 How to Use Getting Started Check Your Skin Tone USING THE SKIN TONE SENSOR Venus Naked Skin® has an integrated Skin Tone Sensor designed to assist you in choosing the appropriate skin tone setting for your skin.

4 1 2 3 5 4 2 The treatment area should be clean, hair-free and dry. 1. Remove any visible hair before using Venus Naked Skin®. Removal of hair ensures the light reaches the hair follicle during treatment. 2.

Cleanse the treatment area and C pat dry. 3 6 Activating Gel 8 8. User Guide 1 HOW TO USE How It Works The Venus Naked Skin® Intense Pulsed Light (IPL) system is derived from professional IPL technology and helps break the cycle of hair regrowth. The light energy is transferred through the <sup>TM</sup> Activating Gel to target the melanin in dark hair, helping to put the hair to sleep. Skin Tone Sensor 1.

Sensor Button 2. Setting Lights 3. Bulb Life Indicator 4. Sensor Window on the back side **HOW TO USE** To measure your skin: 1. Remove the Skin Tone Sensor from the base unit. The RED Skin Tone Sensor light will go out to indicate it is ready. **BEFORE TREAT RESULTS** Setup To power on the device: 1. Place the device unit on a firm surface. The handset should be securely placed in its cradle on the base unit. D <sup>TM</sup> Activating Gel Venus Naked Skin® comes with <sup>TM</sup> Activating Gel.

This is a clear, water-based gel that enables the light to transfer to the hair during treatment to optimize effectiveness.

[You're reading an excerpt. Click here to read official BRAUN IPL VENUS NAKED SKIN BY BRAUN user guide](http://yourpdfguides.com/dref/5484263)  
<http://yourpdfguides.com/dref/5484263>

Do not use Venus Naked Skin® if you do not have Activating Gel. The gel is available for purchase at your retailer. 2. Plug the power supply into an electrical outlet. The Skin Tone Sensor button will illuminate RED to indicate STANDBY mode. A skin tone reading is required to begin. 2. Place the sensor window of the Skin Tone Sensor on the clean, dry area of skin you wish to treat. For comfort and safety, always choose the darkest skin, e.g., for leg, test the outside of the lower leg; for arms, use the back of the hand. For the face, test the cheek or upper lip area. 3. Press the Skin Tone Sensor button to receive a skin tone reading: • If an illuminated WHITE light I (settings I, II, III, IIII) appears, the skin tone in this area is suitable for E treatment.

6. Double check the result by comparing the skin tone of the treated area to the Skin Tone Chart. IMPORTANT: If there is a difference between the skin tone identified by the Skin Tone Sensor and the Skin Tone Chart, you should always use the higher number/darker tone setting (emits less light energy). For example: if the Skin Tone Sensor indicates (III) and the Skin Tone Chart indicates (IIII), you should use (IIII). 7.

Press the Skin Tone Sensor button P repeatedly until the correct skin tone setting appears. STANDBY Mode When the Skin Tone Sensor is placed in the base unit, after 10 minutes of inactivity, the handset will revert to STANDBY mode (indicated by a RED light). A new skin tone reading is required to re-activate the unit. IMPORTANT: Only use the Activating Gel specifically designed for use with this device. In order to prevent gel from drying, larger areas should be treated in sections. 2. Wipe excess gel from your hands. 3. Select a spot on the edge of the treatment area. Place the handset firmly on the skin so that both contact bars on the outer corners of the treatment window are placed in the Activating Gel that is on the skin surface.

A WHITE light signals you have the proper contact. If the white light does not appear, be sure to adjust the handset so that both bars remain in contact with the skin. CORRECT • If an illuminated RED light appears, this indicates the skin tone in this area is not suitable for treatment. I – IIII on the Skin Tone Sensor. HOW TO USE Treatment Before your first treatment, we advise you to perform a patch test on a small area of skin on a body part you plan to treat. Please refer to page 37 for Patch Test details. Once you have established the correct skin tone setting and prepared the treatment area, you are ready to begin your treatment. 1. Apply a thick layer of Activating Gel A to the area. DO NOT RUB IN.

The gel should remain in a thick layer on the skin to ensure proper contact with the Handset. HOW TO USE 4. IMPORTANT: Before treating another body part, you must perform a new measurement. NOTE: If all lights illuminate, this indicates that the measurement taken is not valid, e.g., Skin Tone Sensor G does not have good contact. • Setting I is the HIGHEST power S setting and is for people with fair skin as shown on the Skin Tone Chart—Skin Tone 1. • Setting II is the MEDIUM power S setting and is for people with medium skin tone as shown on the Skin Tone Chart—Skin Tone 2. • Setting III is a LOWER power S setting and is for people with darker skin as shown on the skin tone chart—Skin Tone 3. • Setting IIII is the LOWEST power setting and is for people with darker skin as shown on the Skin Tone Chart—Skin Tone 4.

NOTE: You can choose a higher number/darker tone setting (emits less light energy) than your Skin Tone Sensor measurement result, but a lower number/lighter skin tone setting cannot be selected. Gel layer Skin INCORRECT Gel layer Skin 4. Press the Activation button once. Look away from the handset to avoid the bright light. You will hear a click and may feel a warm sensation. 5. Replace the Skin Tone Sensor in the base unit. The Skin Tone Sensor button will illuminate WHITE which indicates that the unit is ready for use. The device setting will automatically adjust to the Skin Tone Sensor reading. The white setting light indicates the power setting.

3. Remove the handset from your skin. R The two contact bars and treatment window will leave marks in the Activating Gel to help track your treatment progress. 6. Venus Naked Skin® should be used on the skin in continuous rows. Your next selected treatment spot should be adjacent to the previous. There should be no gaps between treatment areas. When placing the handset onto the gel on your skin, the edge of the treatment window and the two contact bars should line up with the marks from the previous treatment area. Facial Treatment (women only) M Using Venus Naked Skin® to remove unwanted facial hair requires care and special guidance to prevent any possible adverse reaction. It is recommended that before treating the face, other body areas are treated first.

Avoid tanning including outdoor sun bathing, tanning treatments in a salon before and during treatment, or use of a tanning lotion or cream before and during treatment as the device may not work for you. Further, darker skin may absorb too much light energy, which may cause discomfort/pain.

IMPORTANT: • Do not treat the eyebrows or areas around the eyes, the scalp, the nostrils, the ears. • Do not use on moles or any dark (pigmented) areas. Before Treatment 1.

Determine if the face may be treated, by using the skin color chart and the skin tone sensor. Because the face is unavoidably exposed to sunlight, it may be darker than other body areas. 2. Perform a patch test: test your skin for reaction in a hair-free area below the ear (see chapter titled “Patch Test”). 3. After this patch test, wait 48 hours, examine the treated area and if it appears normal, begin treating the upper lip, chin or jaw. Getting started 1. Remove all cosmetics, lotions and creams (including self tanning) from the area you wish to treat. 2. Remove any visible hair. 8. When you have finished treating an area, wipe away any excess Activating Gel from the handset and the area you have treated. Place the handset back on the base unit. 9. Check the treated area to see if you develop any skin reactions.

If you experience a side effect, refer to page 29 in this user guide. SLEEP Mode Timeout • When the unit is left “inactive” (i.e., no buttons are pushed or the Skin Tone Sensor is not removed/ replaced) for 55 minutes, then the unit will enter a low power SLEEP mode. • In SLEEP mode, all the lights are off except for the RED Skin Tone Sensor button which will make a single flash once every 3 seconds. You must take a skin tone reading to get out of SLEEP mode. 3. Use the skin tone sensor to confirm suitability for treatment. Do NOT use if the skin tone sensor indicates red (IIII). Always start with the lowest setting when using on the face, regardless of the skin tone sensor reading.

[You're reading an excerpt. Click here to read official BRAUN IPL VENUS NAKED SKIN BY BRAUN user guide](http://yourpdfguides.com/dref/5484263)  
<http://yourpdfguides.com/dref/5484263>